

SIMCRO™

ACCURUS™ INJECTOR

SELF-TENTING

The Accurus self-tenting injector has a radically different design. It provides greater safety and accuracy, saving money and time.



SIMCRO™

ACCURUS™ INJECTOR

SELF-TENTING



For an animal treatment to be effective it needs to be put in the right place. That is where the Accurus out performs any other injector.

The Accurus family of injectors gives the user the ability to perform one-handed, subcutaneous injections. Ensures operator safety while eliminating undue pressure on the animal to activate.

NEEDLE GUARD

The protective needle guard ensures the risk of self-injection is minimised even when the Sekurus is not in use.

REDUCED RISK OF NEEDLE BREAKAGE

A major issue with existing injectors is the needle being the first point of contact with the animal. This results in animal stress and often violent reactions. Because the Accurus needle guard provides the first point of contact (whilst creating the required subcutaneous skin tent) the animal is reassured and settled at the point of injection. This eliminates the commonplace problem of broken or bent needles that create a significant hazard. An added advantage of the needle guard is the protection it gives the needle from surface contamination.



INCREASED OPERATOR SPEED

Administration speed is greatly increased by proficient operators.

Because the Accurus is so simple to use, even untrained staff can easily and quickly achieve correct drug placement.



ONE-HANDED SUBCUTANEOUS INJECTION


Administering a one-handed subcutaneous injection is far easier and less hazardous for the operator than two-handed techniques, particularly in confined spaces such as a cattle crush. Accurus was developed to ensure accurate single-handed subcutaneous drug placement. This allows the user to maintain a safer distance from the animal.



SUBCUTANEOUS INJECTION

The key to efficacy is precision. This is where Accurus is unrivalled. Simcro's exclusive, patented self-tenting needle guard ensures accurate subcutaneous drug placement.

Improved drug placement accuracy = improved drug efficacy.

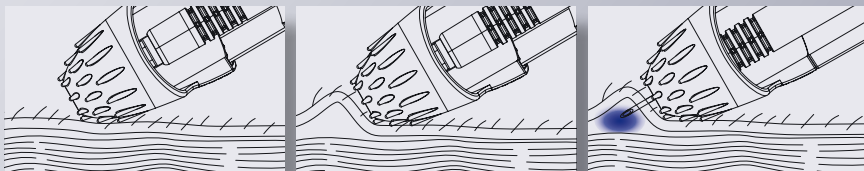

AVAILABLE IN 12, 5, 2 & 0.5 ML SIZES OF TUBE AND BOTTLEMOUNT VARIANTS

The Accurus range includes needle guards designed for perpendicular (intramuscular) and subcutaneous injection. Dose variants include multi-stop selectable 1mL, 2mL, 5mL, 6mL and 12mL options with other dose variants under development.



PARTNER.

INNOVATE.

DELIVER.

Simcro Limited
6 Latham Court
Hamilton 3204
New Zealand
TEL +64 7 846 7723
FAX +64 7 846 7738
EMAIL deliver@simcro.com
www.simcro.com

